



TwinPrep 5™
MANUAL GRINDER/POLISHER

*POWERFUL, RELIABLE & DURABLE
DUAL PLATEN POLISHER FOR
MANUAL SAMPLE PREPARATION*



The **TwinPrep 5™** grinding/polishing machine is designed for manual sample preparation. The powerful, quiet motor drives both platens simultaneously and provides constant high-torque output throughout the speed range, in either 10-500 or 20-1,000 RPM configurations. The 7" color LCD touchscreen with keypad entry is used to control all functions, with an intuitive interface optimized for efficiency. Five (5) user-defined presets allow quick and easy RPM selections, in either clockwise or counterclockwise rotation.

Durable: Engineered and constructed using precision machined, hard-anodized aluminum and stainless steel components to provide maximum longevity.

Easy to use: Straightforward controls and user-friendly functions for quick, repeatable operation.

Either standard or magnetic 8"/203 mm or 10"/254 mm platens may be used with any plain/adhesive backed or magnetic system disc.

Standard/magnetic platens, accessories and consumables are sold separately.

Item	Description
5-7650	TwinPrep 5™, 100-240 V, 10-500 RPM
5-7650-M	TwinPrep 5™, 100-240 V, 20-1,000 RPM

Dims: 31" W x 26" D x 9.5" H (780 x 660 x 240 mm)
Weight: 136 lb. (62 kg)

Features:

- Variable platen speed: either 10-500 RPM or 20-1,000 RPM (10 RPM increments) with five (5) adjustable presets
- 7" LCD touchscreen with keypad entry to control all functions
- Intuitive interface optimized for productivity and function
- Powerful, quiet 0.5 HP (375 W) DC motor that provides constant high-torque output throughout the speed range
- Single motor driving both platens simultaneously
- Clockwise/counterclockwise platen rotation
- Sturdy RIM, aluminum and stainless steel construction
- Molded one-piece corrosion and impact resistant enclosure
- Overflow diversion to protect internal components
- Quick-change platen design, anodized to resist wear and corrosion
- Low profile, tabletop design
- Electronic coolant control with adjustable flow control valves
- Retractable coolant nozzle that allows quick and easy sample/bowl cleaning
- CE compliant for EU
- Two (2) year warranty
- Designed and manufactured by Allied in the USA



Retractable Nozzle



Retractable coolant nozzle allows quick and easy sample/bowl cleaning.

LCD Touchscreen



The 7" color LCD touchscreen is used to control all functions and is extremely easy to navigate, allowing greater efficiency among users.

Standard System

Standard Platens



Standard aluminum platens are precision lapped parallel, then hard anodized for maximum wear resistance. They are used with adhesive backed discs, double sided adhesive discs and paper holding bands, or when converting to a magnetic platen using magnetic bases.

Item	Description
5-2005	8" (203 mm)
5-6010	10" (254 mm)

Disc Holding Bands



Disc holding bands are made of anodized aluminum and secure plain back discs to a standard platen. A protrusion around the top of the band allows variable sized discs (± 3 mm) to be secured.

Item	Description
125-10011	8" (203 mm)
125-10012	10" (254 mm)

Double Sided Adhesive Discs (Pk/10)



Double sided adhesive discs are adhered to a standard platen and used to secure plain back silicon carbide abrasive discs. Each provides 25-75 disc pulls, depending on grit size, cycle time and force.

Item	Description
52-10005	8" (203 mm)
52-10003	10" (254 mm)

Magnetic System

Magnetic Platens



Magnetic platens combine a standard aluminum platen with a cross-hatched magnetic surface for maximum flatness and adherence. Only one platen is needed throughout the grinding and polishing process, allowing all magnetic system cloths/abrasive discs or ferromagnetic discs to be utilized.

Item	Description
5-2005M	8" (203 mm)
5-6010M	10" (254 mm)


Ferromagnetic Support Discs



These reusable ferromagnetic stainless steel discs feature a nonslip backing and allow adhesive backed cloths/abrasive discs to be changed quickly and easily when using magnetic platens.

Item	Description
55-51000	8" (203 mm)
55-51002	10" (254 mm)

PTFE Ferromagnetic Support Discs

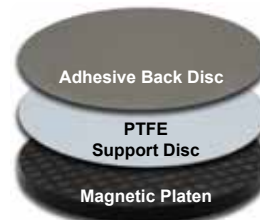

These reusable, durable discs combine a PTFE surface with a ferromagnetic stainless steel substrate. This allows adhesive backed abrasive discs/polishing cloths to be easily removed after use without any adhesive transfer. The nonslip backing keeps the disc in place even for coarse grinding applications.

Item	Description
55-52000	8" (203 mm)
55-52002	10" (254 mm)

Global Distribution of the World's Finest Products



Brochure #: TP5-2020

Equipment is subject to change due to innovations/improvements at any time.



Designed and Manufactured by
Allied High Tech Products, Inc.
in the USA

ALLIED
HIGH TECH PRODUCTS, INC.

2376 E. Pacifica Place, Rancho Dominguez, CA 90220
310.635.2466 Worldwide / 800.675.1118 U.S. & Canada
Visit us online! AlliedHighTech.com